



DeltaCut Biopsy system

For the extraction of biopsy material from organs.

Flexible in use, **quiet** during application.
Simplest handling – sure, precise and fast.


Handy PAJUNK biopsy system “DeltaCut”, into which the biopsy cannula may be inserted in the locked or the unlocked position. Inserting is carried out quite simple, fast, and above all, safe.





Therefore, the perfect positioning of the biopsy cannula is ensured at any time. Exact supervision is possible by means of the 2-colored functioning-status indicator.


The DeltaCut biopsy system will guarantee a faultless and fast extraction of biopsy material.

Optimal extraction of biopsy material through precise biopsy cannulas from PAJUNK.

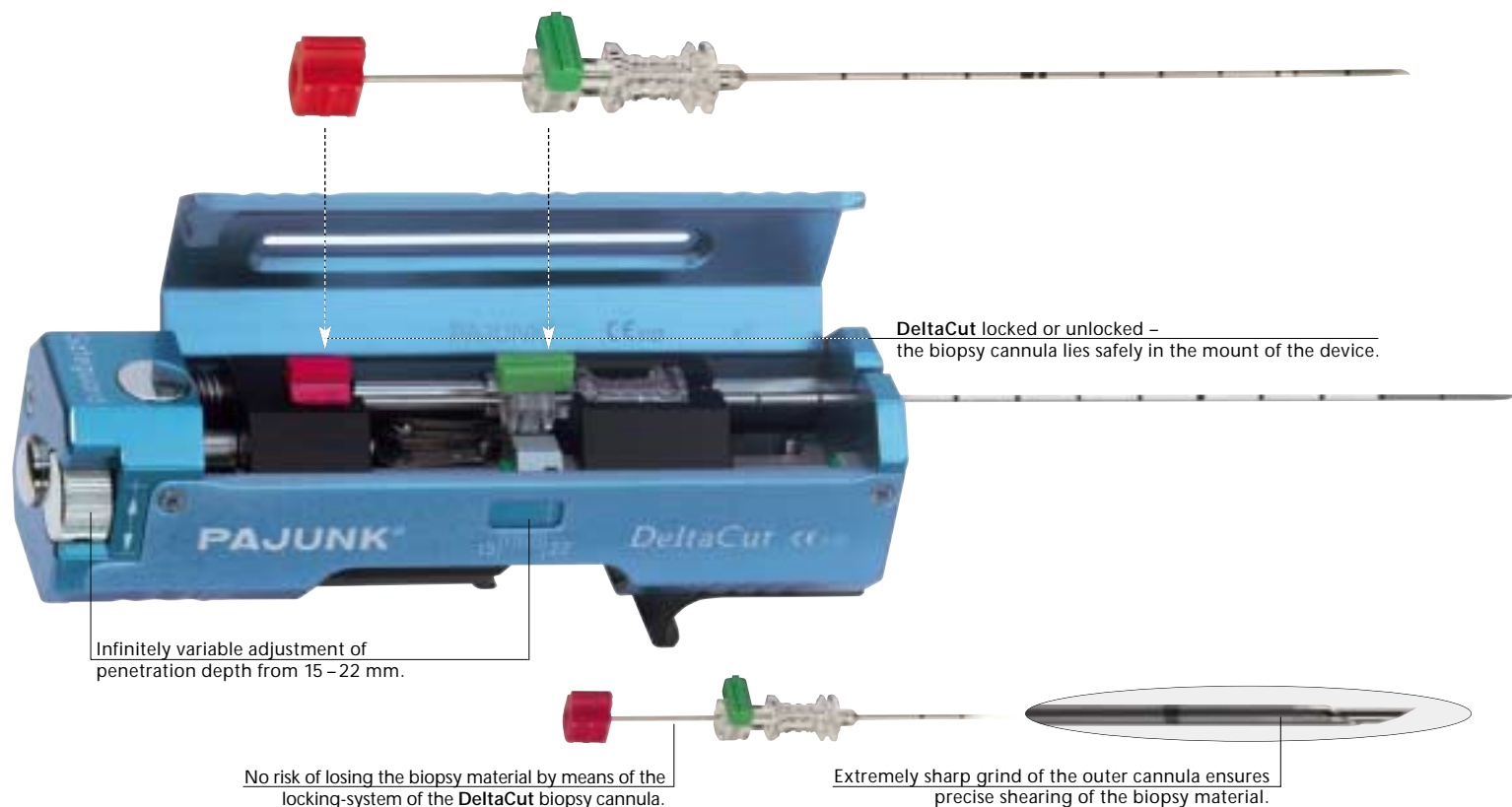Biopsy cannula with the **DeltaCut** grind, exact, precise, with a large biopsy chamber and special locking system for the absolutely safe extraction of biopsy material.

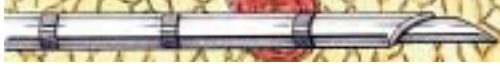

1. Puncture  Biopsy cannula with the **DeltaCut** grind, displaying an exactly sharpened cannula tip. Optimal manufacture: electro-polished and graduated.
2. Triggering  Very large biopsy chamber.
3. Protecting  Automatic closing of the biopsy chamber. Loss of the biopsy material is therefore impossible.
4. Extraction  Cross-section of the biopsy material. Sample will remain in the original spatial condition.

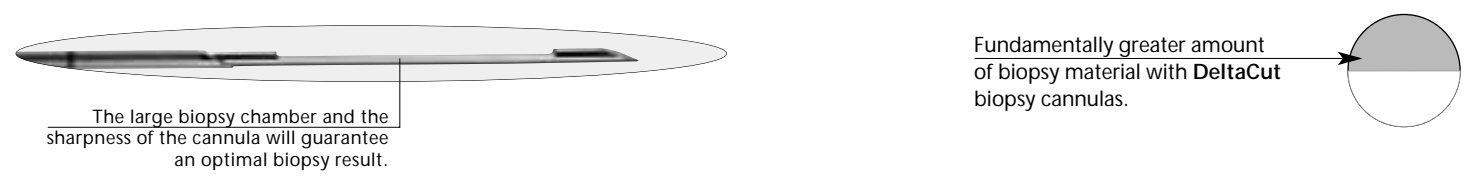


PAJUNK®

DeltaCut Biopsy system

Safe and precise biopsy

- Infinitely variable adjustment of penetration depth from 15–22 mm
- Functioning-status indicator for system supervision
- 2 alternative trigger buttons for ergonomic handling
- Single and multiple biopsies are possible
- Contamination-safe enclosure of the biopsy material within the DeltaCut-cannula

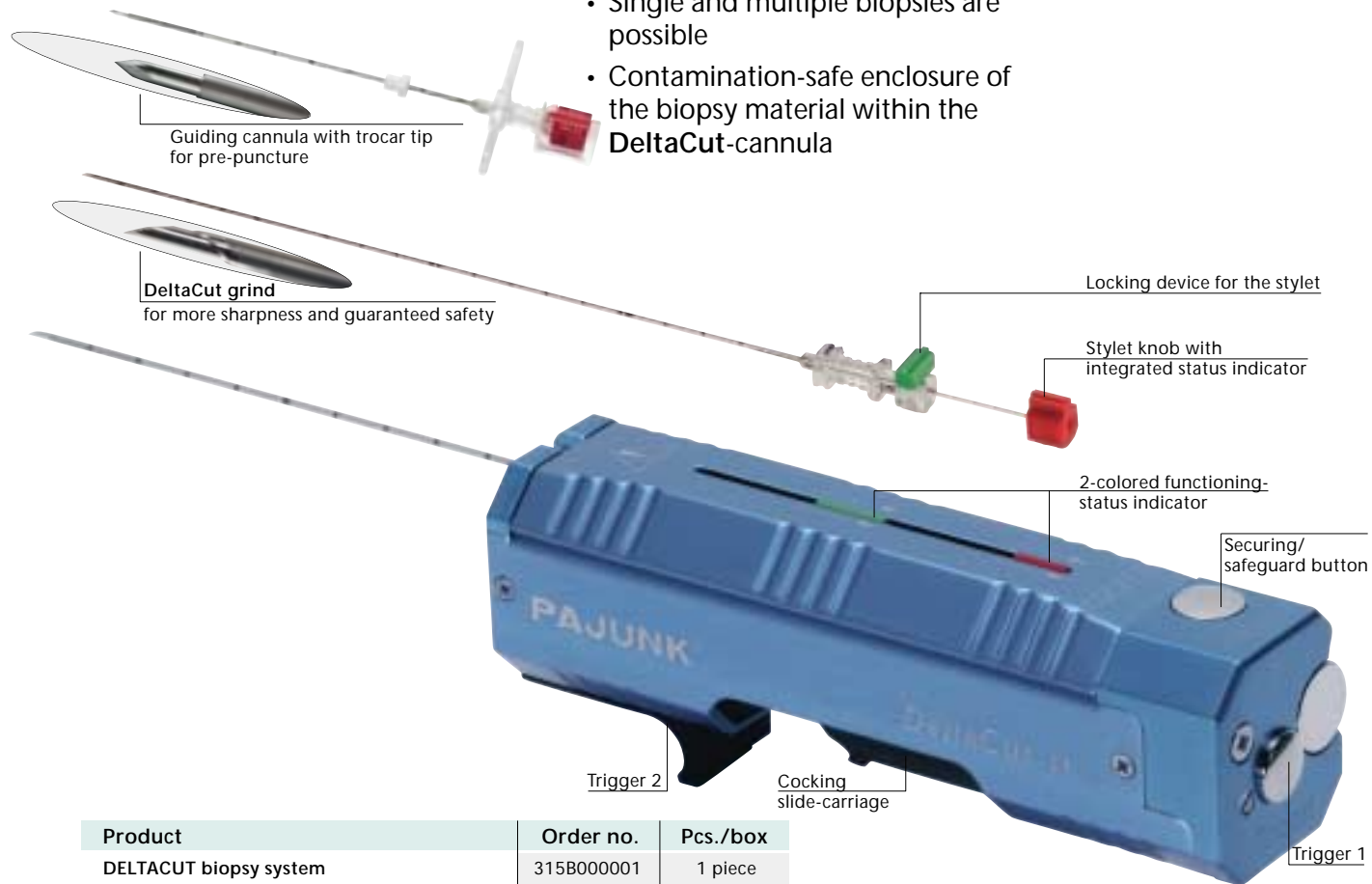

Product	Order no.	Pcs./box
DELTACUT biopsy system	315B000001	1 piece
DeltaCut biopsy cannula, 18G (1.20) x 100 mm	315S120100	10 pieces
DeltaCut biopsy cannula, 18G (1.20) x 150 mm	315S120150	10 pieces
DeltaCut biopsy cannula, 18G (1.20) x 200 mm	315S120200	10 pieces
DeltaCut biopsy cannula, 18G (1.20) x 280 mm	315S120280	10 pieces
DeltaCut biopsy cannula, 16G (1.60) x 150 mm	315S160150	10 pieces

Product	Order no.	Pcs./box
Guiding cannula with trocar tip		
Ø 1,50 x 44 for biopsy cannula Ø 1.20 x 100	313S150014	10 pieces
Ø 1,50 x 94 for biopsy cannula Ø 1.20 x 150	313S150064	10 pieces
Ø 1,50 x 144 for biopsy cannula Ø 1.20 x 200	313S150114	10 pieces
Ø 1,50 x 224 for biopsy cannula Ø 1.20 x 280	313S150194	10 pieces
Ø 2,00 x 94 for biopsy cannula Ø 1.60 x 150	313S200064	10 pieces

Organ-biopsy

PAJUNK®

DeltaCut Biopsy system

*Safe, quick and
precise biopsy*



Organ-biopsy